

Dispermill X-Proof HX400

The Dispermill X-Proof HX400 is designed for the bigger laboratory research and development projects and pilot plants. The H-shaped floor stand ensures stability on the dissolver and the counterweight provides a balance that eases to move the machine up and down with the ergonomic lever. The powerful motor maintains its selected speed, independent of changes in viscosity. This guarantees optimum performance for your dissolving process. Besides, the stainless steel design is easy to maintain.

SAFETY FIRST

Your safety is very important to us. Safety interlocks for vessel detection and operation height are developed to ensure that no hazardous situations occur. This means that you can focus on your activities without being worried. The X-Proof HX400 is designed to operate in potentially explosive atmospheres.

APPLICATION

The Dispermill X-Proof HX-400 fills the gap between the laboratory and the production dissolvers. Do you need a bigger testing sample or a smaller production batch and a machine that operates in potentially explosive atmospheres? That's where the Dispermill X-Proof HX400 comes in.

FEATURES

- ✓ Speed control by frequency inverter
- ✓ Digital speed readout
- ✓ Separate electronic cabinet
- ✓ Vessel clamping system
- 5 m cable between dissolver and electric cabinet included (other length on request)
- 150 mm dispersing disc included





Dispermill X-Proof HX400


Order code

4.0 kW, 0 - 4.000 RPM, 15 L - 150 L

10-1102

N21027 www.dispermill.com



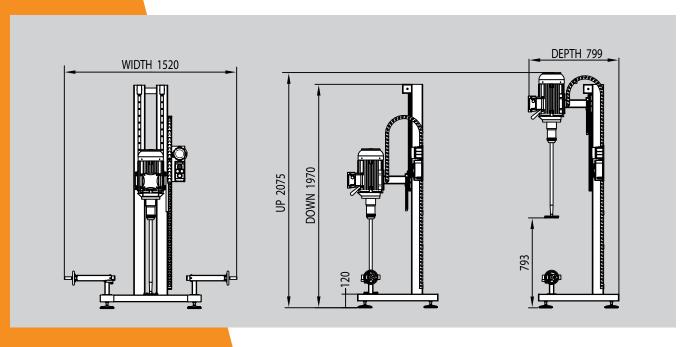

TECHNICAL SPECIFICATIONS


Power	4,0kW / 5 HP
Supply	400 V +/- 10% 50/60 Hz, 16 Ampere
Protection class	(Ex) II 2G IIB T3
Speed	0 till 4.000 RPM
Thermal protection	PTC Temperature sensors
Volume range vessel	15 - 150 L
Diameter vessel	150 - 630 mm (min max.) ¹
Height vessel	300 - 780 mm (min max.)
Dispersing disc	150 mm (interchangeable)
Dimensions (in mm)	799 x 1520 x 1970/2075 (D x W x H)
Weight	+/- 170 KG
Readout	True shaft speed on clear blue display
Safety CE	Conformity Directive 2006/42/EC
Electrical cabinet	Dimensions: 500 x 400 x 210 mm. Requires placement outside of ATEX zone.

¹ in combination with default clamping syster

OPTIONS





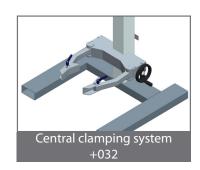

+015 Power readout in kW

+011 Mechanical timer with automatic stop 120 mir

+050 Vessel cover systems (on request)

www.dispermill.com